

Gluteus Maximus Injection Simulator



Size

12.5 x 13.6 x 7.8 (in)

Set Includes

1 main body

- 1 stand for Hochstetter method
- 1 controller
- 1 set of sensor pads

This simulator is designed to practice in the three ways of injection into the buttocks. It is made of a special resin that allows this model to be repeatedly injected without leaving a trace. The equipped indicator shows the accuracy of needle tip positioning and depth.

Features


- 1. The model is made from three individual layers : which are skin, muscular and skeleton. It allows identifying the injection spot.
- 2. The model is designed for three methods: dorsogluteal site, Hochstetter's point and Clark's point.
- 3. Attached indicator shows the accuracy of injection. If the injection spot is correct, the syringe can inject liquid.
- 4. It allows palpating the landmarks such as greater trochanter, superior anterior iliac spine, posterior superior iliac spine and sacrum.
- 5. Both skin and the muscle are replaceable.

Training Skills

· Gluteus maximus injection

Replacement Parts

11213-010

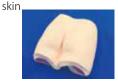

11213-020 sens<u>or for quarter-trich</u>otomy (pair)



11213-030 sensor for hochstetter method (pair)



11213-040 sensor for clark method (pair)



11213-050 controller

