




**Puerperal lower torso with anatomically correct landmarks and elastic abdominal wall facilitates realistic training in assessment of uterine involution in early puerperium.**

**Features**

1. Anatomically correct landmarks for palpation including the pubic symphysis, the umbilicus, ribs and the uterus
2. Soft and elastic abdominal wall that permits highly realistic palpation to assess the degree of uterine involution and to measure the fundal height
3. Four variation of uterine inserts that are interchangeable

**Training Skills**

- Assessment of iterus in early puerperium
- Fundal height measurement by a measure tape
- Palpation of uterine fundus
- Perineal cleansing and assessment
- Demonstration of postpartum massage
- Assessment of anal prolapse

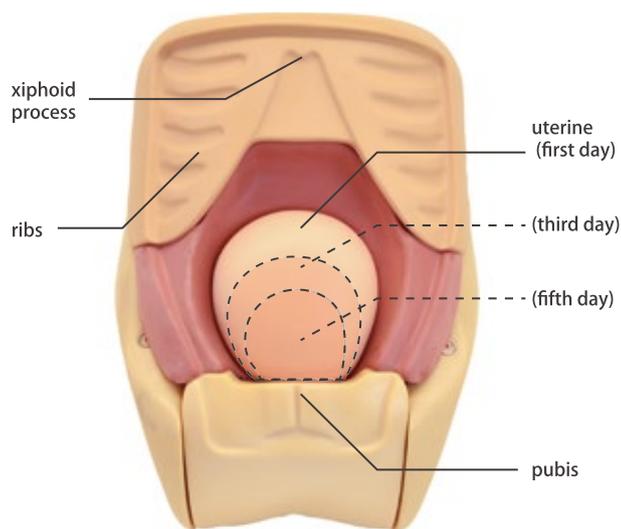

**Variation of uterine (four)**

first day normal	first day abnormal	third day normal	fifth day normal
---------------------	-----------------------	---------------------	---------------------

**Normal** (boundary between abdominal wall and uterine is hard and clear)      **Abnormal** (uterine is soft or unclear)

**Palpation of uterine fundus**



**Measurement of fundal height**



**Specifications :**

Set Includes:		Manikin Size:	Replacement Parts:
1 puerperal lower torso	1 talcum powder	W28 x D42 x H26 cm	(11418-010) 1 skin cover (for uterine involution assessment)
1 set of uterine module (4 variations)	1 instruction manual		(11418-020) 1 set of uterine module (4 variations)